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INSTRUCTIONS MANUAL  
INFRARED THERMOMETER  
FOR OBJECT'S SURFACE  
MODEL FH04

TECHNICAL SPECIFICATIONS

Model : FH04  
The rated power: 2XAAA 1,5V =  
Power : When off 10uV  
When measurement 30 mW  
Operation Distance : 5 ~ 15 Mts

PLEASE READ THE INSTRUCTIONS CAREFULLY OF THIS MANUAL BEFORE  
STARTING.

DO NOT HANDLE THIS INFRARED THERMOMETER INTERNALLY.


The texts, photos, colors, figures and data correspond to the current technical level  
from the moment they have been printed. We reserve the right to modifications  
motivated by the permanent development of the technique in our

Product introduction and classification  
 This infrared thermometer is a thermometer that measures the temperature of the human body using the principle receiving an infrared signal. When used, only forehead detection is required to measure the body temperature quickly and accurately.  
 Product Classification:   
 Basic principle  
 Any object with a temperature above absolute zero will radiate a certain proportion of energy infrared according to its own temperature. The quantity of the radiator energy and its distribution by length of wave are closely related to the object surface temperature. According to this principle, the temperature of the forehead surface can be measured accurately, and the body temperature exact can be determined and displayed by compensating for the difference between the surface temperature of the forehead and actual body temperature.

- Product characteristics
1. Specially designed to measure body temperature, with dynamic compensation for room temperature and body temperature.
  2. Exclusive use of infrared probe for temperature measurement, with high measurement precision and more stable performance.
  3. Provides the sound notification function of higher body temperature.
  4. Provides the sound notification function of higher body temperature.
  5. Capable of storing 32 sets of measurement data. 6.3-color backlit LCD digital display (red, orange, green)
  7. Providing two temperature modes, ie Fahrenheit scale and Celsius scale, for selection.
  8. Provides automatic shutdown function to save electric power.
  9. Small size, reasonable structure and convenient operation.
  10. Select body temperature measurement mode to measure temperature of any object with a temperature below 100 ° C and its emissivity equal to 0.95

**THERMOMETER PARTS:**

1. Battery cover
2. IR sensor
3. LCD screen
4. Body / Object Mode or less
5. C / F or Add mode
6. Memory query or advanced configuration
7. Measurement trigger
8. Grip
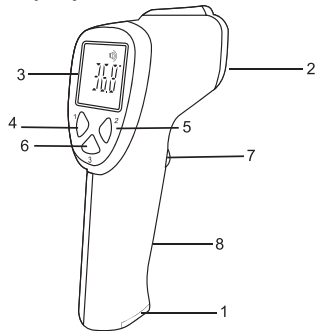
| Function definition | Icon                                                                                | Details            |                                                                                                                    |
|---------------------|-------------------------------------------------------------------------------------|--------------------|--------------------------------------------------------------------------------------------------------------------|
| Level batteries     |  | When is Visible    | The battery is in the low lever but the the thermometer works correctly<br>Replace the battery as soon as possible |
|                     |                                                                                     | When Blinks        | The battery is exhausted and the thermometer cannot work correctly. Replace the battery immediately                |
|                     |                                                                                     | When is Invisible  | The battery power is sufficient                                                                                    |
| Laser               |  | When is Visible    | Laser on (when the machine installed a laser head)                                                                 |
|                     |                                                                                     | When is Invisible  | Laser off (when the machine installed a laser head)                                                                |
| Measurement mode    |                                                                                     | Body Temperature   | Body Mode                                                                                                          |
|                     |                                                                                     | Object Temperature | Object mode                                                                                                        |

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0197

RoHS

|                 |      |                              |
|-----------------|------|------------------------------|
| Reading scale   | °C   | Celsius reading              |
|                 | °F   | Fahrenheit reading           |
| Monitor reading | 8888 | Temperature value            |
| Memory          | M    | Memory value                 |
| Buzzer          |      | Buzzer On                    |
|                 |      | Buzzer Off                   |
| Bluetooth       |      | Bluetooth on / of (optional) |

DEFINITION OF THE FUNCTION OF ALL BUTTONS.

| No. Buttons |            | Description                                                                                                                                                                                                         |
|-------------|------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| 1           |            | Multifunction button, short press to open / close backlight, press it for 3 seconds to change the body mode to object mode or object mode to body mode                                                              |
| 2           |            | Multi-function button, short press to turn on or turn off the buzzer                                                                                                                                                |
| 3           | Set memory | Multifunction button, use with 1 and 2 button, press it briefly to check the memory data, press the test button to exit, press it for 3 seconds, enter advanced settings, use with button 1 and button 2 as part 7. |

SETTINGS

| Menu | Function                   | "+"            | "-"            | Default | Comments                  |
|------|----------------------------|----------------|----------------|---------|---------------------------|
| F1   | Value adjustment deviation | Decrease 0,1°C | Increase 0,1°C | 0,0°C   | Object mode it's useless. |
| F2   | Value adjustment alarm     | Decrease 0,1°C | Increase 0,1°C | 38°C    | Effective range 3°C       |
| F3   | Change C/F                 | °C             | °F             | °C      |                           |

NOTICE

- The temperature in the human body mode is obtained from the dynamic compensation of the ambient temperature and the surface temperature of the forehead.
- The object temperature mode is to test the surface temperature of an object. The temperature obtained from the forehead in this mode is the temperature obtained from the forehead in this mode the value is simply the temperature of the forehead surface, but not the body temperature.
- The deviation setting can adjust the measurement value from -3.0 °C to 3.0°C according to the test distance, ambient temperature, skin difference, etc. The default value is 0.0 ° C

MEASUREMENT

Body temperature

- Press the measurement button to turn on the thermometer and the startup screen will appear (fig. 9.1). After POST and two beeps, room temperature value will be displayed and you are ready for measurement (fig. 9.2).
  - Make sure the thermometer is in body mode.
  - Maintain a distance of 50mm to 150mm from the upper eyebrows to the probe (Fig. 9.3).
- Press the measurement button and when it beeps the measurement is complete and the value will be displayed (fig. 9.4). If the measurement value exceeds the alarm value (Default value is 38 ° C) it will give "bi. Bi. Bi" as an indication.
- After measurement, if the thermometer is inactive in 20 seconds

GUARANTEE DOCUMENT

Thank you for purchasing this COMELEC product, this device has been checked in all its components meticulously, so we guarantee the correct operation, but for this to be the case we need that YOU read the instruction book carefully, as this will contribute to the best performance of this article.

P.G.E, S.L guarantees to the buyer of this product, a total guarantee of TWO YEARS, governed by the established conditions whose regulation is found in Royal Legislative Decree 1/2007 of 16 December, which approves the revised text of the General Law for the Defense of Consumers andUsers, (articles 114 to 127)

P.G.E, S.L is committed to repairs are completely free, labor and spare parts, including.

The Technical Assistance Service will resolve any repairs that may be required by this apparatus.

Carefully review the INSTRUCTION BOOK included, on the use, maintenance and conservation of the device.

Any anomaly that could occur due to the non-observance of such instructions will not be protected for the guarantee.

The date of sale of the device, for the purposes of warranty period, will be that which appears on the purchase invoice.

Guarantee operation: As regards all those compliance defects detected by the technical assistance, which are manifested within six months of the delivery of the goods, it is assumed that they already existed on that date, the intervention will be completely free, (unless this hypothesis is incompatible with the nature or nature of the conformity defect).

During the following eighteen months, it must be the user who demonstrates. That the lack of conformity It was coming at the time of delivery.

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Post Code, 46970 City: Alaquas VALENCIA

We declare under our responsibility that the product:

Product Description: Infrared thermometer

Model ; FH04

Brand: COMELEC

It complies with the European Union legislation in compliance with the 2004/108 / EU guidelines

The standards apply:

ICC60601-1

IEC60601-1-2

ISO13485

ISO149/1 ISO80601-2-56

Name: Pablo Gonzalez

Position : Administrador

Date: 10 Junio 2020

Signature:

INFORMATION ON RECYCLING

Disposal of the old appliance, the European Directive 2002/96 / EC on waste equipment Electrical and Electronic (WEEE), establishes that old appliances can not be thrownas municipal waste without normal classification. Appliances have to be collected separately with theto optimize the recovery and recycling of the materials they contain and reduce the impact on human health and the environment.

The symbol "container with wheels" crossed out on the product reminds you of obligation, that when you dispose of the device must be collected separately. Consumers should contact the local or retail authority to obtain information about the correct disposal of your old device.

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DATE OF SALE

Seller's signature and stamp

IMPORTANT: To make the guarantee effective, the SELLER must fill out, sign and stamp this Certificate. To receive assistance under the conditions set forth in this certificate, please submit to the Technical Service Authorized, this document next to the Purchase Invoice.

MEASUREMENT OF FIGURES

(fig9.1 )


(fig9. 2)

(fig9. 4)

(fig9. 5)

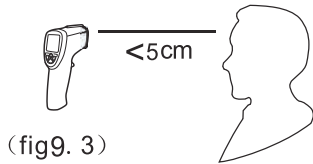(fig9. 3)

 **DANGER**

1. Keep the sensor and probe cavity clean before and after use.
2. To ensure the accuracy of the measurement, it is recommended to start the measurement after ten minutes when transporting the thermometer to a new environment.
3. Wait 10 minutes to measure body temperature after measuring temperature objects extremely high or extremely low temperature.
4. Wait 5 minutes to start a measurement when a target measurement (human objector) is from a environment with huge temperature difference
5. Breeze, water, sweating, cosmetics on the forehead may affect measurement. Do not measure body temperature in 30 minutes after exercise, bath or meals.
6. Do not allow young children to operate the appliance or use it as a toy.

**OBJECT TEMPERATURE**

- Press the measurement button to turn on the thermometer (fig. 9.6).
- After measuring the ambient temperature, make sure that the thermometer is in object mode. (fig.9.7) • (fig.9.8)
- Maintain a vertical distance of 50mm to 150mm from the object to the probe. Press the button and when you give it a "bi ..." measurement finished and the value will be displayed. After measurement, if the thermometer is inactive in 20 seconds

(fig.9.6)

(fig.9.7)

(fig. 9.8)

(fig.9.9)

 **NOTICE**

- 1.The value in this mode is the temperature of the object's surface instead of the core temperature.
2. The default infrared emissivity is 0.95. The reading will deviate from the actual temperature according with different emissivity For example, the reading on stainless steel will obviously be lower than the actual temperature. BE CAREFUL FOR ESCALD.

**EXCEED THE MEASURING RANGE**

Body mode:

When the measurement value is less than 32.0 ° C, the display shows Lo (Fig 9.10) and gives "beep.beep.beep.beep".

When the measurement value is above 42.9 ° C, the display shows Hi (fig 9.11) and gives "bi.bi.bi.bi".

Object mode:

When the measurement value is less than 0 ° C, it displays Lo (fig 9.12) and gives "beep.beep.beep.beep".

When the measurement value is greater than 100.0 ° C, it displays HI (fig 9.13) and gives "bi.bi.bi.bi".

 **AVISO**

When the ambient temperature is below 10°C or above 40°C. Shows ERR. It is not allowed to measure or accuracy is not assured ...

(fig9.10)

(fig9.11)

(fig9.12)

(fig9.13)


(fig9.14)

**BATTERY CHANGE**

Open the battery cover and take out the exhausted battery.

Put 2 AA batteries and close the battery cover. After installing the new battery, the thermometer will emit "Neep. Beep" If there are no beeps, check if the positive and negative poles are correct (see Fig 9.1)

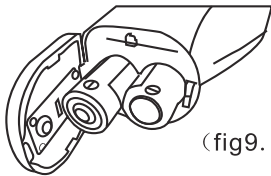(fig9. 1)

 **NOTICE**

Remove the battery in case the thermometer is not used for a long period. Do not expose fire the battery.

Dispose of the battery in accordance with local regulations.

**MAINTENANCE AND TIPS**

Make sure the sensor and capacity are clean, otherwise it will affect the accuracy

Probe cleaning method.

Use a cotton swab or a soft cloth with water or alcohol to clean the case.

Use a cotton swab or a soft cloth with water or alcohol to gently clean the surface sensor or probe cavity. Do not use the thermometer before vaporizing the alcohol.

Please read this manual carefully before use. Make sure the battery is properly installed.

It is not allowed to place the thermometer in any liquid or exposed to strong sunlight or temperatures extremely low.

A strong knock or blow to the product will cause its damage.

Do not dismantle this thermometer yourself.

Keep the thermometer out of the reach of children

Do not use the thermometer under circumstances of strong electromagnetic interference.

Measurement results are likely to fluctuate due to improper forms of measurement.

Practice appropriate measures to improve your skill.

Measurement results cannot replace a doctor's diagnosis.

Special maintenance is unnecessary for this thermometer.

Contact the dealer or manufacturer in the event of a malfunction.

**PROBLEM SOLVING**

| Description                                       | Solutions                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                 |
|---------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| LCD screen "LO" or "HI"                           | <ol style="list-style-type: none"><li>1. Breeze, water, sweating, cosmetics on the forehead may affect the measurement.</li><li>2. Check the offset value adjustment. The default value is 0.8°C.</li><li>3. Whereas if the ambient temperature test changes enormously or if the thermometer is used directly from high temperature injection up to very low temperature, measurement difference will occur. The thermometer should be kept in a relatively stable environment for 10 minutes to obtain balance before starting a new measurement.</li><li>4. Make sure the measurement distance is 5cm to 8cm</li></ol> |
| There's no answer pressing the button measurement | <ol style="list-style-type: none"><li>1. Remove and reassemble the battery.</li><li>2. Check if the thermometer is in the menu settings on the menu setting procedure, the thermometer cannot measure and, therefore it does not respond</li></ol>                                                                                                                                                                                                                                                                                                                                                                        |
| Without display or wrong screen                   | Take out the battery and install the battery again                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                        |
| Without buzzer                                    | Check if the buzzer is off                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                |
| Shut down just after turning on                   | Check the battery level or remove and install it again                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                    |