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GUARANTEE DOCUMENT

Thank you for purchasing this COMELEC product, this device has been checked in all its components meticulously, so we guarantee the correct operation, but for this to be the case we need that YOU read the instruction book carefully, as this will contribute to the best performance of this article.

P.G.E, S.L guarantees to the buyer of this product, a total guarantee of TWO YEARS, governed by the established conditions whose regulation is found in Royal Legislative Decree 1/2007 of 16 December, which approves the revised text of the General Law for the Defense of Consumers andUsers, (articles 114 to 127)

P.G.E, S.L is committed to repairs are completely free, labor and spare parts, including.

The Technical Assistance Service will resolve any repairs that may be required by this apparatus.

Carefully review the INSTRUCTION BOOK included, on the use, maintenance and conservation of the device.

Any anomaly that could occur due to the non-observance of such instructions will not be protected for the guarantee.

The date of sale of the device, for the purposes of warranty period, will be that which appears on the purchase invoice.

Guarantee operation: As regards all those compliance defects detected by the technical assistance, which are manifested within six months of the delivery of the goods, it is assumed that they already existed on that date, the intervention will be completely free, (unless this hypothesis is incompatible with the nature or nature of the conformity defect).

During the following eighteen months, it must be the user who demonstrates. That the lack of conformity It was coming at the time of delivery.

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Post Code, 46970 City: Alaquas VALENCIA

We declare under our responsibility that the product:

Product Description: Infrared body thermometer

Model ; NDT 0

Brand: COMELEC

It complies with the European Union legislation in compliance with the 2004/108 / EU guidelines

The standards apply:

EN 55014-1: 2006+A1:+A2:

EN 55014-2: 1997+A1:+A2:

EN 61000-3-2:2006+A1:+A2:

EN 61000-3-3:2013

Name: Pablo Gonzalez

Position : Administrador

Date: 10 Enero 2.019

Signature: 

INFORMATION ON RECYCLING

Disposal of the old appliance, the European Directive 2002/96 / EC on waste equipment Electrical and Electronic (WEEE), establishes that old appliances can not be thrownas municipal waste without normal classification. Appliances have to be collected separately with theto optimize the recovery and recycling of the materials they contain and reduce the impact on human health and the environment.

The symbol "container with wheels" crossed out on the product reminds you of obligation, that when you dispose of the device must be collected separately. Consumers should contact the local or retail authority to obtain information about the correct disposal of your old device.

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DATE OF SALE

Seller's signature and stamp

IMPORTANT: To make the guarantee effective, the SELLER must fill out, sign and stamp this Certificate. To receive assistance under the conditions set forth in this certificate, please submit to the Technical Service Authorized, this document next to the Purchase Invoice.

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Table2-Guidance and MANUFACTURER'S declaration -electromagnetic IMMUNITY-for all EQUIPMENT and SYSTEMS

Guidance and MANUFACTURER'S declaration -ELECTROMAGNETIC IMMUNITY			
The [EQUIPMENT or SYSTEM] is intended for use in the electromagnetic environment specified below. The customer or the user of the [EQUIPMENT or SYSTEM] should assure that it is used in such an environment.			
IMMUNITY test	IEC 60601 Test level	Compliance Level	Electromagnetic environment guidance
Electrostatic discharge (ESD) IEC 61000-4-2	± 6 kV Contact ± 8 kV Air		Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30 %.
Electrical fast transient / burst IEC 61000-44	± 2 kV for power supply lines ± 1 kV for input/output lines		Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	1 kV) line(s) to line(s) g ± 2kV line(s) to earth		Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	< 5% U _T (> 95% dip in U _T) for 0,5 cycle (40% U _T (60% dip in U _T) for 5 cycles 70% U _T (30% dip in U _T for 25 cycles (> 95% dip in U _T) for 5 s		Mains power quality should be that of a typical commercial or hospital environment. If the user of the [EQUIPMENT or SYSTEM] requires continued operation during power mains interruptions, it is recommended that the [EQUIPMENT or SYSTEM] be powered from an uninterruptible power supply or a battery.
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	3A/m		Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.
NOTE U _T is the a.c. mains voltage prior to application of the test level.			

10

19

INSTRUCTIONS MANUAL INFRARED BODY THERMOMETER MODEL NDT 0

SUITABLE FOR FAMILY USE.
This product is for home use only.

TECHNICAL SPECIFICATIONS

Model : NDT 0
The rated power : 2XAAA 1,5V =
Power : 1,5mW
Operation Distance : 3 ~ 5 Mts

PLEASE READ THE INSTRUCTIONS CAREFULLY OF THIS MANUAL BEFORE STARTING.

DO NOT HANDLE THIS INFRARED BODY THERMOMETER INTERNALLY.

The texts, photos, colors, figures and data correspond to the current technical level from the moment they have been printed. We reserve the right to modifications motivated by the permanent development of the technique in our

ROW	Guidance and MANUFACTURER'S declaration -ELECTROMAGNETIC EMISSIONS EMISSIONS: for all EQUIPMENT and SYSTEMS		
1	Guidance and declaration of the MANUFACTURER - ELECTROMAGNETIC EMISSIONS		
2	The [EQUIPMENT or SYSTEM] is intended for use in the electromagnetic environment specified below. The customer or the user of the [EQUIPMENT or SYSTEM] should assure that it is used in such an environment.		
3	Emissions test	Compliance	Electromagnetic environment - guidance
4	RF Emissions CISPR 11	Group 1	The [EQUIPMENT or SYSTEM] uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
5	RF Emissions CISPR 11	Group 2	The [EQUIPMENT or SYSTEM] must emit electromagnetic energy in order to perform its intended function. Nearby electronic equipment may be affected.
6	RF Emissions CISPR 11	Class (A or B)	
7	Harmonic emissions IEC 61000-3-2	[Class A, B, C, D or Not applicable]	
8	Voltage fluctuations/ flicker emissions IEC 61000-3-3	[Complies or Not applicable]	
9		[See 6.8.3.201 a) 3)and Figure 201]	The [EQUIPMENT or SYSTEM] is suitable for use in all establishments, including domestic establishments and those directly connected to the public lowvoltage power supply network that supplies buildings used for domestic purposes.
10		[See 6.8.3.201 a) 3)and Figure 201]	The [EQUIPMENT or SYSTEM] is suitable for use in all establishments other than domestic, and may be used in domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes, provided the following warning is heeded: Warning: This equipment/system is intended for use by healthcare professionals only. This equipment/ system may cause radio interference or may disrupt the operation of nearby equipment. It may be necessary to take mitigation measures, such as re-orienting or relocating the EQUIPMENT or SYSTEM] or shielding the location.
11		[See 6.8.3.201 a) 3)and Figure 201]	The EQUIPMENT or SYSTEM] is suitable for use in all establishments other than domestic and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
12	RF emissions CISPR 14-1	Complies	The [EQUIPMENT] is not suitable for interconnection with other equipment.
13	RF emissions CISPR 15	Complies	The [EQUIPMENT] is not suitable for interconnection with other equipment.

INTRODUCTION

Thank you for purchasing the infrared body thermometer. Infrared Thermometer TI7150 is mainly used for patient's body temperature measurement. Usually it is operated by medical staffs, patient himself or his families. And it cannot be used for life support or rescue. It will not cause any accident even if the device is invalid. It cannot be used for special or emergency measurements. Suitable for all different ages, infrared body thermometer is used for monitoring human body temperature. This infrared thermometer can be also used to measure the temperature of objects such as bath water, milk or others. Its non-contact design makes this infrared thermometer safer and cleaner than most of other thermometers. This infrared thermometer can even be used while a child is sleeping.

'Intensive measurement technology' has been adopted, which is with exclusive patent. The thermometer will measure many times one second. The measuring data will be more correct because the average value will be displayed. Please read the following instructions carefully before operation. With proper care and use, it will be reliable for many years.

SAFETY PRECAUTIONS

In order to enable you use the thermometer safely and properly and keep away from the harm because of misusing, please note the warning signs and symbols in the manual. Warning signs, symbols and their meanings are explained as follows:

	Warning	Casualties may occurred when misuse
	Note	Misuse may cause personal injury or product damaged.
		Operating instructions

SYMBOLS

	Person/Objetc/ModE		Memory
	Scan		Keep Dry
	Frágile, Handle with care		This way up
	Fabricante		Standby
	Authorised Representative In The Comunity European		
	European Union' s Waste Electrical and Electric Equipment in accordance with Article 11 (2) Of Directive 2002/961EC WEEE) To protect the environment, product must be disposed at the definite places according to your national or local regulations		

SPECIFICATIONS

Measuring Range:

Body: 34,0°C ~ 42,1°C (93,2°F ~107,7°F)

Object: 5°C ~100°C (41°F ~ 212°F)

Accuracy:

Body 0,2°C (36,0°C ~ 39,0°C)
0,3°C (34,0°C ~36,0°C o 39,0°C ~ 42,1°C)
Objeto: 1°C

Measuring position: Forehead

Memory: 10 readings

Operation distance: 3 ~ 5Cm

Operating environment : 5°C ~ 35,0°C 85%RH
700HPA ~ 1060HPA

Transport/Storage 25°C ~ 55,0°C 95%RH

environment 500HPA ~ 1060HPA

Dimension: 134x48x31 mm

Expected service life: 2 Years

Power: 2xAAA 1.5V alkaline batteries

Power consumption 1. 5mW

Battery Life 1 Year

Weight: 76g (without battery)

WARRANTY CONDITIONS:

- 1) Repair the thermometer for free within one year from the date of purchase if used correctly. Service: The thermometer should be returned to manufacturers and repaired by professionally trained people.
2. In the warranty period, if the units fail, you can return it to the store you bought or to our headquarters along with this Product Warranty.
3. The following situation cannot be repaired for free:
 - Products damaged by misuse or user review without written factory permission.
 - Products damaged by collision or fall.
 - Products damaged by force majeure, such as fires, earthquakes, floods, and other abnormal voltages.
 - Products damaged by failure do not follow instructions.
 - Remove or remove the no series and repair labels from the products.
 - Consumables
 - No failure of the product itself
4. Take care of this Product Warranty, no additional warranty card will be issued if lost.
5. We repair the thermometer free of charge in accordance with this Product Warranty instruction. However, this warranty does not restrict the legal rights of customers.

Note:

The manufacturer will provide circuit diagrams, descriptions of component parts lists, Calibration instructions to assist the designated service person to repair the parts.

COMELEC PRODUCT WARRANTY

- 1) Free repair within one year from the date of purchase if used correctly.
2. In the warranty service request, you must have a stamp of the date of purchase and the place of purchase (including name and address). Be sure to buy this product with the store seal on this Warranty.
3. Out of warranty, repair service will be charged.
 - Do not reissue, please keep this warranty safely.


Taking a Reading - Object

- 1. If the thermometer is off, press the 'Power' button to switch the thermometer on.
- 2. Press button to set the mode of object measurement. The screen should appear the follow picture.
- 3. Aim the front of the thermometer at the object. Hold the thermometer approximately 3 to 5 cm.
- 4. Short press the button. The screen should appear the follow picture:
- 5.Heading beep , the temperature measurement value will be displaded for 2 seconds. The display will illuminate for 3 seconds.
- 6.The thermometer will automatically beep and shut-off after approximately one minute of non-use.


Memory Recall.

- 1. Press the ' ' button to recall up to 10 measurement results stored in memory.

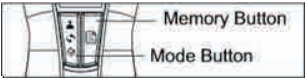- 2. Thea symbol or symbol will indicate a person or object temperature which was taken.
- 3. The happy face or sad face symbol will appear 5 seconds after viewing the measurement results stored in memory.

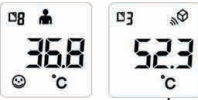
NOTE:

After 10 measurement results have been stored in the memory, the oldest reading will be automatically deleted.

Clearing Memory

- 1. Press the Power button to turn the thermometer on.
- 1. To clear the memory, simply hold the button down until the thermometer beeps and two horizontal lines appear.

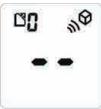Cleaning & Storage

To ensure accurate measurements, clean the lens before and after each use. Fingerprints and dirt can cause inaccurate readings.

To Clean

- Use an alcohol swab or a soft cotton swab moistened with alcohol to clean the thermometer lens.
- Do not use paper towels, because the fibers may scratch the lens.
- Do not use water to clean the lens or thermometer.
- Only use a soft dry cloth to clean the thermometer
- Make sure that liquid won't enter the interior of the thermometer.
- Never use abrasive cleaning agents, thinners or gasoline to clean.
- Do not immerse into water or other cleaning liquids.

TO STORE

- Remove the battery if the thermometer will not be used for an extended period of time.
- Store in a sun-free and dust-free dry area at a temperature between -25 °C —55 °C (-13 °F ~ 131°F) with relative humidity not exceeding 95%.
- Protect the thermometer from severe impact or shock. Keep stored in the protective storage bag.

WARNING

The measurement result of this IR thermometer is not intended as a substitute for consultation with your physician.

WARNING

It is very dangerous if the patients make self-judgment and treatment according to the measurement results. So please get the treatment according to the doctor's guidance.

Self-judgement may worsen the illness. The thermometer is only for forehead temperature. So don't use for other parts (such as in the ear).

Otherwise, the thermometer may hurt the body.

If the infrared probe is damaged or dirty, it may cause incorrect measurement results.

When the infrared probe is dirty, please wipe with soft and dry fabric gently.

Keep the thermometer in measure location for 30 minutes before operation if there is a big difference in temperature between storage site and measuring site.

Do not measure body temperature within 30 minutes after meal, sports or bath.

Otherwise measurement results may get incorrect. Please note the battery polarity when replace the battery

Otherwise it may cause leakage, heat, rupture even damage to the thermometer.

The thermometer should be used in working temperature range, please check the product label.

Otherwise the infrared thermometer may not work or make the measurement incorrect. The operator cannot touch battery and patients at the same time

Please pay attention to the angle of the thermometer's front part to prevent hands from scratching

IMPORTANT NOTE

Keep the infrared body thermometer far from the place where kid can touch.


The measurement result cannot be substituted for the consultation from doctors

Keep product in place away from direct sunshine, high temperature, heavy humidity, dust, fire-side, easy shock and vibration areas.

Please do not use the thermometer in electromagnetic interference environment

UNIT

OVERVIEW

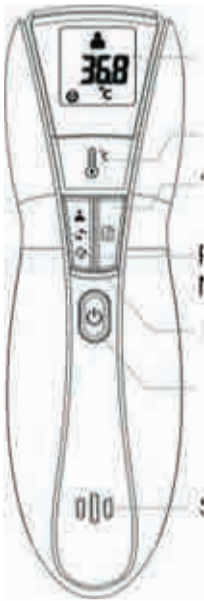LCD Display

SCAN Button


MEM Button
(memory)

Person/Object
Mode Button

Led light

Power Button

Speaker

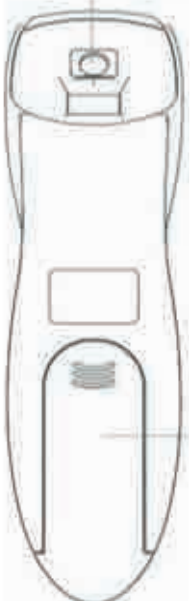Infrared Lens

Cover Battery

FRONT


REAR

LCD Display	MEANING	DESCRIPTION
°C/°F	Fahrenheit Celsius	
	Forehead Scan Mode	
	Object scan mode	
	Memory symbol and measurement mode	
	Scanning in progress	
	Low Battery	Unit will work;replace batteries soon
	Happy Face	Indicates temperature is lower than 104°F.38°C
	Sad Face	Indicates temperature is over 104°F.38°C
	Measurement is below normal temperature range.	Retake temperature ensuring proper positioning
	Measurement is above normal temperature range	Retake temperature ensuring proper positioning
	Operating temperature is outside of the permissible range (32°-104°F/ 0-40°C)	Take thermomter to a room where the temperatuer is within the permissible range. Wait 45 minutes to stabilize thermometer .
	Blank display	Check battery polarity.May need to replace batteries.
	Low battery warning	Replace the batteries immediately

REPLACING THE BATTERY

When the battery symbol appears on the screen, the batteries will need to be replaced soon. When the letters 'LO' and the battery symbol appears on the display, the batteries need to be replaced before taking another measurement. Rechargeable batteries also can be used in this thermometer.

1. Gently slide the battery cover back.
2. Carefully remove the old batteries and properly discard.
3. Insert new batteries according to the proper polarity
4. Slide the battery cover back on.

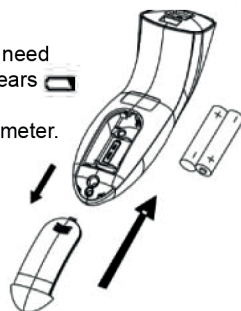NOTE:

Please properly dispose of the batteries away from children and heat.

Please remove the battery if you do not use the infrared thermometer for a long time.

In order to protect the environment, please send the waste batteries to the appropriate collection points according to the national or local laws.

OPERATING GUIDE
Switching Between Person or Object

1. Press the 'Power' button to switch on the thermometer
2. Press and hold both the button to change the measuring mode of person and object.

Object


Body

NOTES:

- a. Anthropometric mode: the measurement results have been amended to the oral temperature.
- b. Objects mode is calibration mode. Calibrate once every two years.

Switching Between Fahrenheit or Celsius

1. Press the 'Power' button to switch on the thermometer.
2. Press and hold both the and buttons for about 3 seconds. This will change the mode to either °F or °C. Once the thermometer beeps, the unit is ready for measurement.

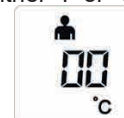Taking a Reading-Person

1. If the thermometer is off, press the 'Power' button to switch the thermometer on.
2. Please remove hair,aim the front of the thermometer upon the person's forehead. Hold the thermometer approximately 3 to 5 cm away from the forehead.

NOTE: Attempting to take a measurement anywhere else on the body may result in an inaccurate reading.

- 3.Short press the button. The screen should appear the follow picture.
- 4.Heading beep ,the temperature measurement value will be displayed for 2 seconds. The display will illuminate for 3 seconds.

NOTE: A happy face symbol will appear on the display if the temperature is below 38.°C (100.4°F). If the temperature reading is above 38.0 °C (100.4°F), a sad face symbol will appear.

If the temperature is outside of the thermometer's sensory range, 'LO' or 'Hi' will appear on the display. However, if the temperature is outside of the thermometer's permissible measurement range, 'Err' will appear on the display.

- 5.The thermometer will automatically beep and shut-off after approximately one minute of non-use.

Helpful Hints

- For accurate results, please do the temperature measurement at least 30 minutes after exercising, bathing or returning from the outdoors.
- Wait at least 1 minute between 2 measurements.
- If multiple readings are taken in a short period of time, use the average of those readings.
- Because everyone's temperature varies slightly, keep a record of your body temperature to serve as a reference when judging a fever.